



Oval, POLYsmooth™



Oval, MESMO®



Round, POLYsmooth™

Gluteus Implants

Silicone-gel-filled gluteus implants available in two shapes and surfaces, with planar or concave back

Implants made by POLYTECH – QUALITY made in Germany



Gluteus Implants, € 0483

Product Description

Silicone gel-filled gluteus implants for compensation of soft tissue deficiencies are available with two different surfaces:

- POLYsmooth™ (smooth)
- MESMO® (micro-textured)

The surgeon can choose among oval and round shapes of various sizes as well as between a concave or planar implant posterior. Thus he is in a position to achieve the optimal result for either augmentation or reconstruction of the gluteus region in patients of both sexes.

Information on Materials

The shell of all gluteus implants consists of a strong, resilient and highly elastic silicone elastomer equipped with a barrier layer. It is filled with a highly cross-linked silicone gel. Due to the high stress that can possibly act on the implants, a gel with firm consistency has been chosen. The surface of the implants is either smooth, POLYsmooth™, or micro-textured, MESMO®. The MESMO® texture features a small pore size of an average 100 to 150µm.

Production of the implants is based on a quality management according to ISO 9001 and ISO 13485. All materials used are of implantable medical quality.

Indications

- Augmentation
- Contour correction
- Reconstruction, e.g. after trauma
- Implant exchange

Important Information

Gluteal implants are delivered individually packed and sterile.


Attention: Please note the detailed instructions for use included with each product.

Gluteus implants – concave posterior


Silicone-gel filled, individually packed, sterile, € 0483

Order n°	Vol. (ml)	Dimensions (mm)			Posterior
		A	B	C	

Oval, POLYsmooth™					
10300-115 ¹	115	80	120	27	concave
10300-215 ¹	215	105	142	30	concave
10300-290 ¹	290	114	152	33	concave
10300-360 ¹	360	125	163	34	concave
10300-450 ¹	450	132	174	38	concave

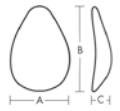

Oval, MESMO®					
15300-115 ¹	115	80	120	27	concave
15300-215 ¹	215	105	142	30	concave
15300-290 ¹	290	114	152	33	concave
15300-360 ¹	360	125	163	34	concave
15300-450 ¹	450	132	174	38	concave

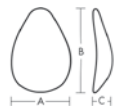

Round, POLYsmooth™					
10301-200 ¹	200	107	107	35	concave
10301-250 ¹	250	112	112	41	concave
10301-300 ¹	300	118	118	44	concave
10301-360 ¹	360	127	127	43	concave




Gluteus implants – planar posterior


Silicone-gel filled, individually packed, sterile, € 0483

Order n°	Vol. (ml)	Dimensions (mm)			Posterior
		A	B	C	

Oval, POLYsmooth™					
10300-240 ¹	240	90	130	40	planar
10300-305 ¹	305	95	139	45	planar
10300-365 ¹	365	100	148	48	planar
10300-420 ¹	420	109	152	50	planar

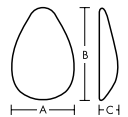

Oval, MESMO®					
15300-240 ¹	240	90	130	40	planar
15300-305 ¹	305	95	139	45	planar
15300-365 ¹	365	100	148	48	planar
15300-420 ¹	420	109	152	50	planar

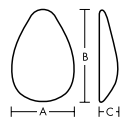

Round, POLYsmooth™					
10301-225 ¹	225	107	107	35	planar
10301-275	275	112	112	41	planar
10301-340	340	118	118	44	planar
10301-390 ¹	390	128	128	44	planar
10301-560 ¹	560	141	141	49	planar



approximate measurements

On stock	¹ Supply after order ≈ 6 weeks
----------	---

